



80
YEARS OF QUALITY
1932 ~ 2012

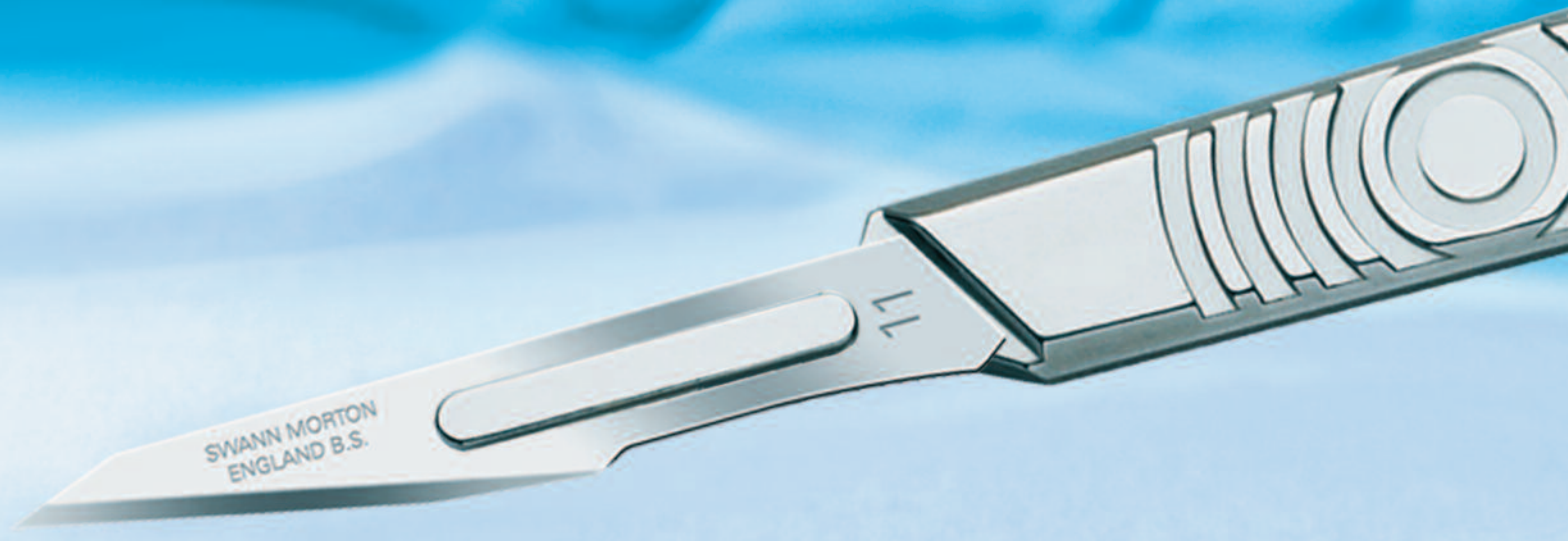

Surgical blades and handles range

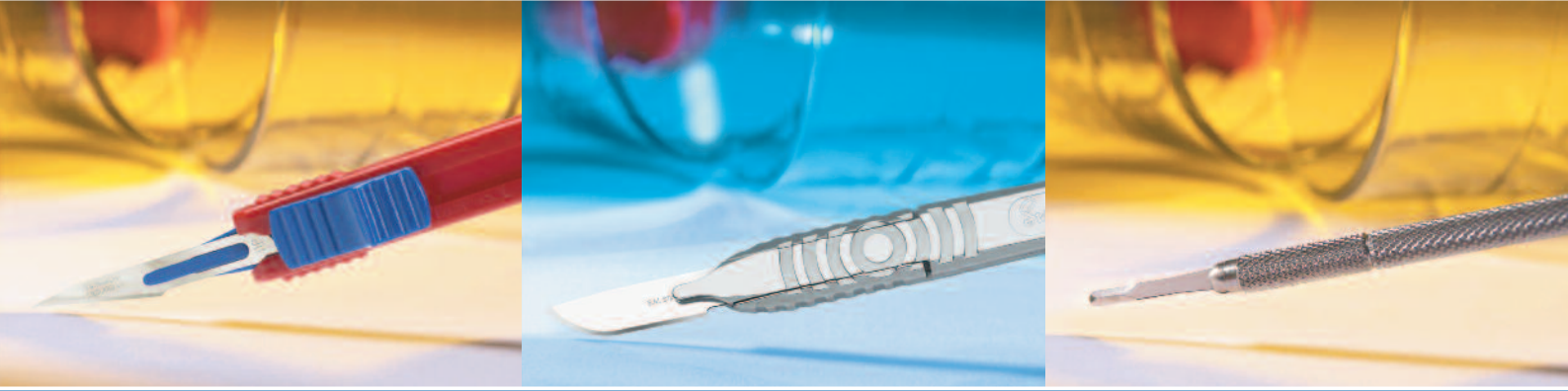

For more information contact:

Swann-Morton®
Made in Sheffield - England

Owlerton Green, Sheffield S6 2BJ
Telephone: 0114 2344231, Sales: 0114 2344223, Fax: 0114 2314966

uksales@swann-morton.com
exportsales@swann-morton.com
orders@swann-morton.com
accounts@swann-morton.com
www.swann-morton.com

The home of quality surgical blades

AWD/02/12/Issue 2

Swann-Morton®
Made in Sheffield - England

The
Quality
Edge



Swann-Morton scalpels and blades have a reputation for quality worldwide that is second-to-none.

That reputation has been built on exceptionally – and consistently – high standards of durability and sharpness. Swann-Morton is accredited to BS EN ISO 13485, enabling the company to CE mark its products under the Medical Device Directive (93/42/EEC). All standard fitting surgical blades and handles are produced in accordance with BS 2982, ISO 7740 / BS EN 27740. Every single batch of blades is subject to the strictest possible quality control procedures at every single stage of its manufacture, and the company also has a fully equipped microbiological laboratory.

Such attention to quality is, of course, vital for Swann-Morton. The company's products play an important role in one of the most demanding environments in the world – the operating theatre.

For the same reason, Swann-Morton blades are subject to constant development and refinement. As surgical techniques become ever more highly refined, so must Swann-Morton develop blades that are equal to the challenge.

Surgeons in just about every country in the world trust other people's lives to Swann-Morton blades. There can surely be no better evidence of the importance the company places on quality.



Scan this code
for more product
information

The Swann-Morton Surgical Blade Remover

Swann-Morton's Surgical Blade Remover is for single use and is available sterile in boxes of 50 and non-sterile in boxes of 100. It is designed to remove surgical blades from standard No. 3 and 4 fitment handles and protects hospital staff removing scalpel blades after operations, guarding against cuts and abrasions.

Skin Graft Blade (SG3)

In stainless steel, with 5cm long cutting edge. Fits a No.3 fitment handle. Sterilised, packed individually, boxed in 20's.

Myringotomy Blade

Long, angled, stainless steel blade with short cutting edges at tip. Purpose-designed for both up and down cutting. Individually packed, sterilised and boxed in 10's. This blade fits the 'FINE' range of handles.

Surgical Blades and Handles by

Swann-Morton®

www.swann-morton.com

Disposable Scalpels

Available with all standard stainless steel blade shapes, complete with CE marked polystyrene handles with a graduated measure (for indication only), individually sealed in a tray designed to assist with safe passing protocols, sterilized and boxed in units of 10.

Guarded Disposable Scalpels

Available with the most popular stainless steel blade shapes, and are fitted with a polyethylene protective guard. The sterile scalpels are individually wrapped in a view pack, and boxed in units of 10. Also available in a bulk, non-sterile format. For use in custom procedure trays and kits.

Sabre Disposable Scalpels

A choice of 3 'Sabre' sterile disposable scalpels – E/11, D/15 and B/23 – each supplied individually packed, boxed in 10's.

Disposable Scalpels with Retractable Blade

These disposables have been specially designed to reduce accident risk potential by having a high-visibility red handle, retractable blade for Safe Passing operated by a colour-coded button which audibly clicks to indicate open and locked position, and ribbed, slip-resistant finger grips for secure hold in gloved hands. They are available with stainless steel blade No's 10 (black button), 11P (blue button) and 15 (green button), in sterile packs of 25.

Also available in non-sterile packs of 500. For use in custom procedure trays and kits.

Standard Surgical Blades and Handles – No.3 Fitment

Surgical Blade No's 6, 9, 10, 10A, 11, 11P, E/11, Sabre E/11, 12, 12D, 13, 14, 15, 15A, 15C, 15T, Sabre D/15, 16, 40 fit No's 3, 3L, 5, 7, 9, B3 & B3L handles.

Handles are available in stainless steel, with selected sizes available in nickel. All handles are individually wrapped. The easily gripped B3 & B3L are light, adaptable stainless handles, supplied individually in plastic wallets.

Standard Surgical Blades and Handles – No.4 Fitment

Surgical Blade No's 18, 19, 20, 21, 22, 22A, 23, 24, 25, 25A, 26, 27, 36 fit No's 4, 4L and 6 handles.

Handles are available in stainless steel, with selected sizes available in nickel. All handles are individually wrapped.

Standard surgical blades in carbon steel (sterile or non-sterile) or in stainless steel (sterile). Boxed in 100's (sterile blades individually foil-wrapped; non-sterile sealed in foil packs of 5). Non-sterile blades also available individually packed in peel-open packets, boxed in 100's. Sabre E/11 and D/15 blades are produced in carbon steel only, sterile, boxed in 100's.

'Fine' Blades and Handles

The 'Fine' range has small, slim, precise blades in stainless steel for delicate operations. Six special handles, including the SF13 and SF23 with independently screwed top sections to facilitate blade-changing, take all blades, supplied in individual metal foil packets, available in dispenser units of 25. Handles are supplied individually in plastic wallets.

Cygnetic® Surgical Blades and Handles

For use within rigorous Orthopaedic and General surgical procedures, Cygnetic® provides the surgeon with additional strength and stability through the patented combination of its specially designed blade and handle.

Handles are available in Stainless Steel, Non-Sterile, individually sleeved. Blades in Stainless Steel, Sterile, individually V.P.I. protected in Metal Foil Packets and boxed in units of 50.

Post-Mortem Knife – PM40 and PM40B

Post-Mortem knife with a choice of heavy-duty carbon steel blades. Screw-type fitting to stainless steel handle. Spare blades, individually wrapped, boxed in 10's.

Post-Mortem PM8 Handle & PM60 and PM60B Blades

Complementing existing Post-Mortem range, offering easy blade change, in non-sterile carbon steel. Both blades V.P.I. protected in single-peel packets, boxed in 50's. Handles supplied individually in plastic wallets.

Stitch Cutters

Requiring no handles, the standard stitch cutter is in carbon steel, the midi and long versions in stainless steel. The No 3 stainless stitch cutter can be used in conjunction with a No 3 handle for more aseptic techniques.

Available individually foil-wrapped in sterile packs of 100.

Sterilisation is by Gamma Radiation (minimum dosage 25kGy, 2.5 Mrad) carried out under strict laboratory control in Swann-Morton (Services) Ltd's Cobalt 60 continuous irradiation plant.

All items are shown actual size



"Swann-Morton" and the 'ring pattern' logo are registered trademarks of Swann-Morton Limited and related companies. The products shown are manufactured for use by trained health care professionals and can cause injury if used or handled incorrectly. The information and illustrations in connection with the products shown on this leaflet are for general guidance only and users must seek proper training and instruction in their use before using or handling such products. Swann-Morton Limited exclude all liability for injury to the extent permitted by law.

Stitch Cutter

Stainless steel with No.3 handle fitment. Individually foil-wrapped in boxes of 100.

For more information about Cygnetic® product refer to our website at: www.swann-morton.com

For cleaning and reprocessing instructions and to view an online video of this procedure, go to <http://www.swann-morton.com/cygnetic>

Cygnetic®
(Showing the
Attach/Detach
position)

Cervical Biopsy Blade

In stainless steel, sterilised, supplied in single plastic trays in boxes of 10. Fits the Long and Short Swann Major Handles which can be viewed on the website along with the Swann-Major blades.